

BIOS living

Electronic Glass Scale

Instruction Manual

SC420



Introduction

Thank you for purchasing the BIOS Living Electronic Glass Scale. This scale was developed with state-of-the-art technology and digital electronics; this precision device has a capacity of 330 lbs / 150 kg, and is equipped with an automatic OFF function.

Features

- 2 $\frac{3}{4}$ " / 7.3 cm LCD readout
- High tempered glass platform
- Battery included (2 x CR2032)
- Displays weight in LBS / KG / ST
- Auto OFF feature
- Division 0.2 lbs / 0.1 kg
- Accurate to 330 lbs / 150 kg

Starting your Scale

1. The lithium batteries (CR2032) are included with your bath scale. Remove the battery cover from the back of the scale, by lifting the tab. Remove the protective film from batteries.
2. To change the unit from LBS / KG/ ST, turn the scale on and press the button to toggle to the desired measurement.

Operation

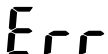
1. Place the scale on a hard, flat surface. Always weigh yourself on the same scale placed on the same floor surface.

NOTE: If the scale is placed on a soft surface, such as a carpet, or on an uneven floor the results may be inaccurate or an error message may occur.

2. Turn ON the scale by tapping the scale in the middle with your foot. When the scale turns ON and 0.0 is displayed on the screen, stand on the scale platform to take your measurement.
3. When the scale flashes your reading 3 times it is finished taking your measurement. The scale will automatically shut off after 10 seconds.

Troubleshooting

Error warnings will appear on the screen, as shown, if the following occur.

Error Displayed	Reason for Error
	Indicates low battery
	If the scale is overloaded - Over capacity (330 lbs / 150 kg / 124 st) remove weight immediately from scale

Advice for Use and Care

- Avoid using platform when wet - platform will be slippery when wet. Keep scale dry.
- Avoid an extremely high or low temperature environment.
- Avoid excessive impact or vibration to the unit.
- Clean the scale with a damp cloth and avoid water getting inside scale. Do not use chemical cleaning agents.
- Do not store the scale in an upright position when not in use as this will drain the battery.
- Remove the battery if the scale is not being used for a long time.

Five Year Warranty

If this product proves to be defective in material or workmanship within five years of purchase, please return it to the address below. It will be repaired or replaced without charge upon receipt of the unit prepaid with \$5.00 to cover handling, packaging and return postage. Please include proof of purchase, your full name, address, daytime phone number or email address.

Thermor Ltd.
16975 Leslie Street,
Newmarket, ON L3Y 9A1
www.biosmedical.com