
Project Type	Date: January 13 2021	Revision: A	Thermor
Instructions	Customer: Thermor Ltd Job Number: TH4720 Item Number: 57006 Designer: Alex Vranjesevic		Contact: Graphic's Department graphics@thermor-ins.com
UPC (For Reference Only) 0 57475 57006 4	Reason For Project: New Item: New Branding: Design Update: X N/A:	Artwork is not to be amended or scaled. If any changes are required please send through Thermor office with instructions, and we will amend here and resend artwork through. **PLEASE RETURN WITH SIGNATURE**	Colour Breakdown K


BIOS living

SITZ BATH
57006
Instruction Manual

BASIN PREPARATION:

1. Raise toilet seat and cover.
2. Place Sitz Bath basin on bowl with "Front" designation towards front of bowl.

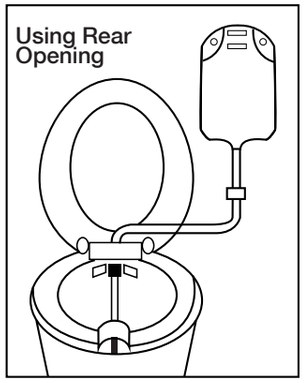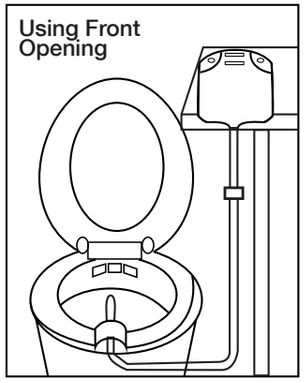3. Fill basin with warm water (38°C to 40°C / 100°F to 105°F) 1/2 to 2/3 full.
4. Insert tubing through front or rear entry hole in Sitz Bath basin. Snap or secure tube into channel in basin bottom.

BAG PREPARATION:

1. Position flow control clamp toward Sitz Bath Bag.
2. Using clamp, close tubing by sliding clamp over tube from wide to narrow side.
3. Open top of bag and fill with hot water (40°C to 43°C / 105°F to 110°F).
4. Place or hang bag above Sitz Bath.

TREATMENT:

1. Sit in Sitz Bath for 15-20 minutes.
2. Maintain water temperature by adding water by opening clamp, (overflow will drain through openings at back of basin).
3. Once desired temperature is reached, close clamp.
4. After treatment, empty and rinse basin.

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM

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