Graphic's Department
graphics@thermor-ins.com

This artwork and design are the property of Thermor Ltd. and are intended solely for the production of authorized products. Reproduction, distribution, or use without written consent is strictly prohibited.

PROJECT DETAILS

Project Type	Customer	Model Number	Designer	Colour Breakdown
Instructions	Thermor Ltd	LF766	Alex Vranjesevic	CMYK
Date Modified	Job Number	UPC - Barcode	Revision	
December 2 2025	TH6859	0 57475 26766 7	D	

DO NOT PRINT (Cut) DO NOT PRINT (Crease)


BIOS living™
Transfer Bench

LF766
User's Manual


Weight Capacity 300 lbs/136kg


1 Check all parts


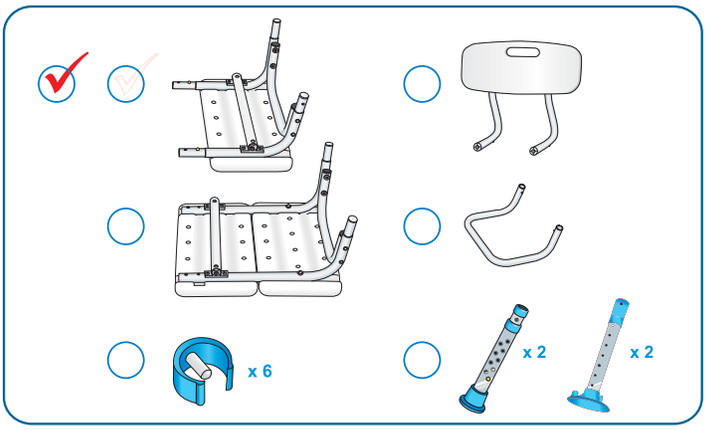
2 Tool-free setup


To make instructions easier to follow, colours are used for contrast and do not reflect the colour of the product.

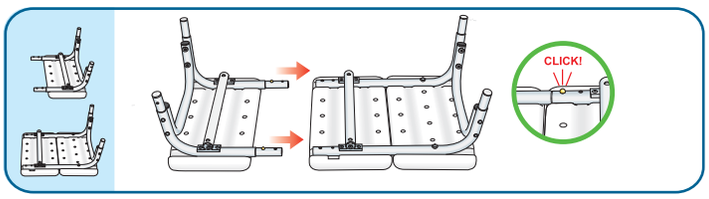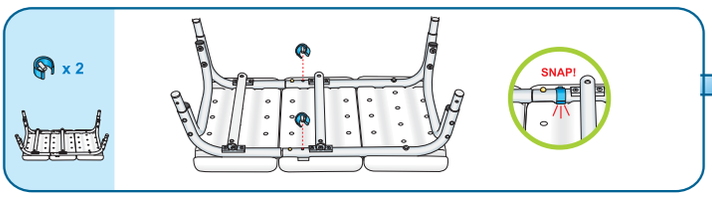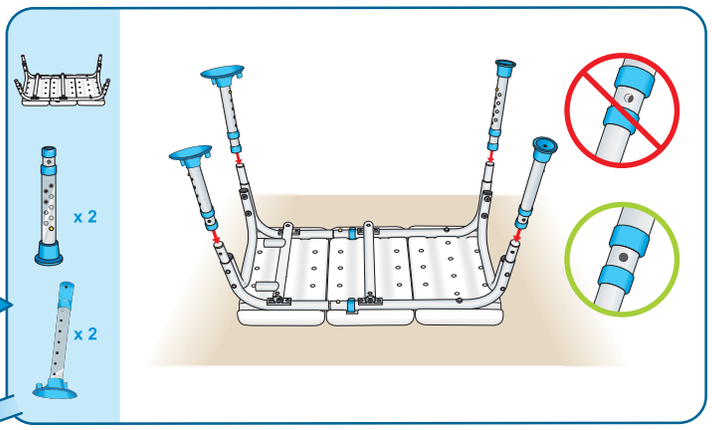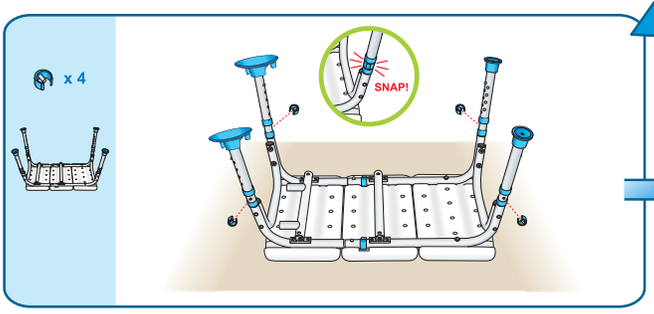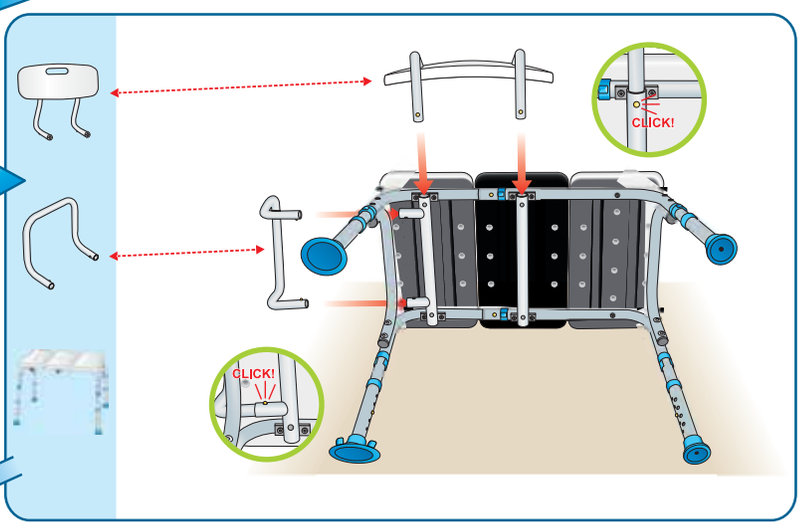MAINTENANCE

1. Use a mild soap and water solution, rinse with clear water, then dry with a soft cloth. Ensure that the bench is well rinsed before using.
2. Never use organic solvents, abrasive cleaners, or scouring pads on any part of the bench.
3. Check fit and tightness of all nuts, bolts, and knobs weekly to ensure that the bench is safe to use.

SAFETY PRECAUTIONS

DO NOT ATTEMPT TO ADJUST OR OPERATE TRANSFER BENCH WITHOUT READING ALL INSTRUCTIONS CAREFULLY.

- Make certain that the height adjustment snap buttons fully protrude through the same respective holes of each leg extension. This will ensure that the leg extensions are securely locked in position and an even height is achieved.
- **DO NOT** stand on transfer bench.
- BIOS Medical assumes no responsibility for any damage or injury caused by improper installation, assembly, or use of this product.

Warranty

Three Year Warranty

This Transfer Bench has a one year warranty to be free of manufacturing defects in materials and workmanship under normal applications for the lifetime of the original owner. If this product becomes inoperable due to defect, return the product with all component pieces and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling, or damage.

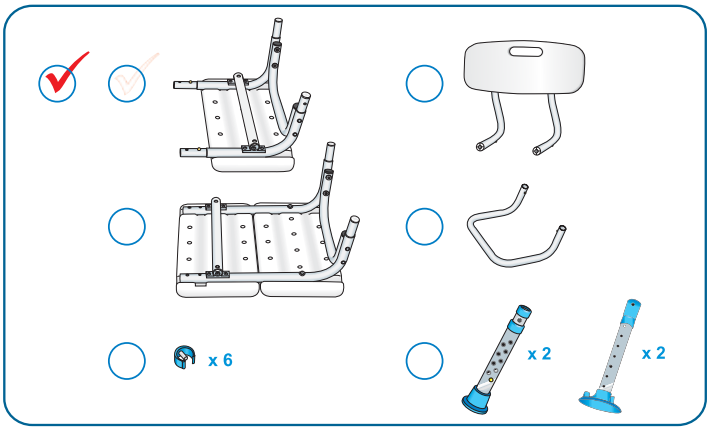
Ship the unit prepaid and insured (at owner's option) to:

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THERMOR LTD.