
Project Type	Date: May 17 2021 Customer: Thermor Job Number: TH6860 Item Number: LF773 Designer: Heather Yoksimovich	Revision: A	Thermor **Contact:** Graphic's Department graphics@thermor-ins.com
UPC (For Reference Only) 0 57475 26773 5	Reason For Project: New Item: **X** New Branding: Design Update: N/A:	Artwork is not to be amended or scaled. If any changes are required please send through Thermor office with instructions, and we will amend here and resend artwork through. **PLEASE RETURN WITH SIGNATURE**	**Colour Breakdown** C M Y K


BIOS living

BATH BOARD ASSEMBLY / ASSEMBLAGE DE LA PLANCHE DE BAIN

1 Place plastic bolts through adjustment slots.

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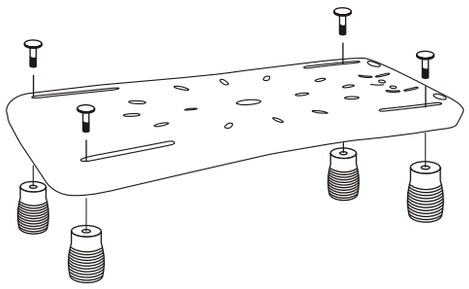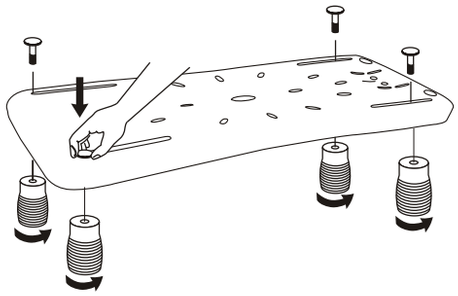
2 Align adjustment leg with plastic bolt and turn counter clockwise until snug. Do not over tighten adjustment legs.

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
1 Place handle into holes on the bath board and align slotted plastic screws.

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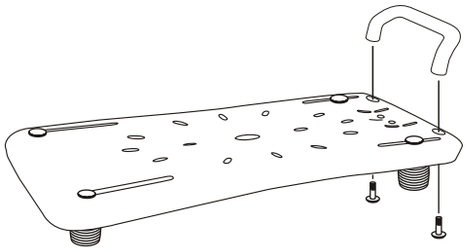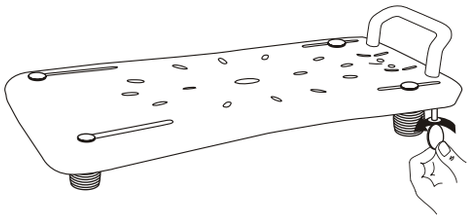2 Tighten slotted plastic screws using a coin until snug. Do not over tighten.

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

1 Loosen adjustment legs slightly and move them outwards so they make contact with the inner walls of the bathtub.

Re-tighten adjustment legs. The bath board should now fit securely in the bathtub. It should not exceed the outside wall of the bathtub by more than 2 cm.

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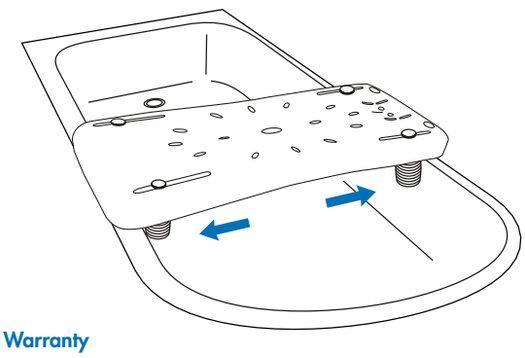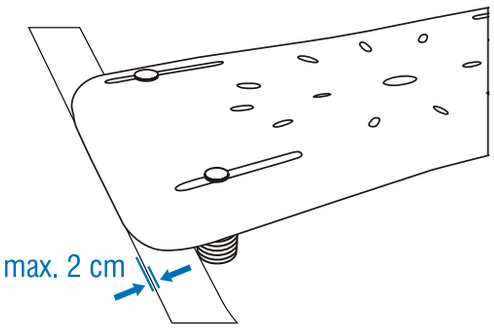One Year Warranty

This Bath Board has a 1 year warranty to be free of manufacturing defects in materials and workmanship under normal applications for 1 year of the original owner. If this product becomes inoperable due to defect and requires repair, return the product with all component pieces and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply if the product is subject to misuse, neglect, rough handling or damage.

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Ship the unit prepaid and insured (at owner's option) to:
Thermor Ltd.
Attn: Repair Department
16975 Leslie Street
Newmarket, ON L3Y 9A1

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