BIOS living

TOILET SAFETY RAIL
LF758
INSTRUCTION MANUAL

- Sturdy steel construction provides stable support
- Adjustable width: 19”- 22” / 48.25 - 55.88 cm
- Weight capacity: 300 lbs / 136 kg

ASSEMBLY INSTRUCTIONS:

Step 1. Assemble armrests by sliding armrests into tubing with corresponding numbers. Use C clips to secure. **(Figure 1)**

Step 2. Slide assembled armrests in bracket that is to be secured to toilet using the push pin holes to adjust width. Make sure the brass pushpins are fully extended and secure. **(Figure 2)**

Step 3. Remove existing toilet seat and save the bolts for reattachment. **(Figure 3)**

Step 4. Place assembled frame on toilet. Position toilet seat on top and secure with bolts. **(Figure 4)**

Step 5. Place the rubber tips under the toilet seat. The weight of the user will keep the Toilet rail in place. **(Figure 5)**

CAUTION: This device is to be used as a safety frame only. It is not designed to support the entire weight of the individual. Periodically check frame screws to make sure rail is securely attached. Tighten if necessary.

THREE YEAR LIMITED WARRANTY

CUSTOMER SERVICE
Toll Free: 1-800-387-8520

Figure 1

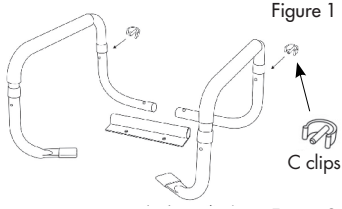


Figure 2

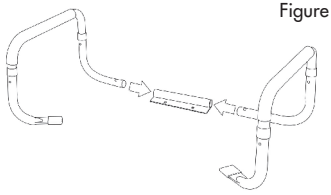


Figure 3

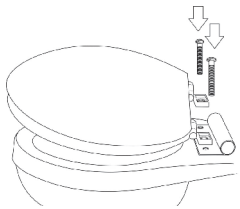


Figure 4


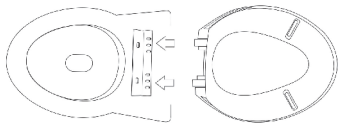


Figure 5

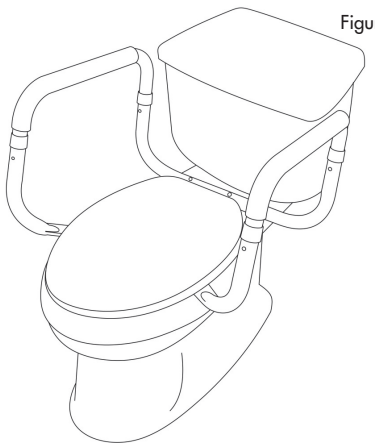

Paper
Papier

BIOS | Medical
MANUFACTURED BY/FABRIQUÉ PAR :
THERMOR LTD.
16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA/FABRIQUÉ EN CHINE
WWW.BIOSMEDICAL.COM
100% Canadian Owned & Operated
Propriétaire-exploitant canadien à 100 %